

MAGNUM[™] EX

2. Generation

Universal Magnet Plate

Powerful Ring Magnets for Demanding Applications

■ Rapid Separation

Uses our strongest grade ring magnets enabling 100% faster bead separation compared to standard magnet plates

■ Large & Small

The Magnum EX magnet handles PCR plates just as well as full volume deep well plates

■ Compatible

Works with full-, semi-, and non-skirted 96-well PCR plates, as well as shallow, half-deep, and deep well plates with U- and V-bottoms.

■ Removable Corner Tabs


Dropped your magnet plate and broke a corner tab? Easy, just pop in a new one!
A corner tab gets in the way of your plate gripper? No problem - remove it and run without!

■ Full Aluminum Construction

Solid metal construction supports workflows where temperature control is desired (up to 80°C)

■ Includes Patented* PAQLOC[®] Spring Lock

In just a few seconds, enable the springs for perfect supernatant removal on the liquid handler, or lock the springs for manual pipetting on the lab bench

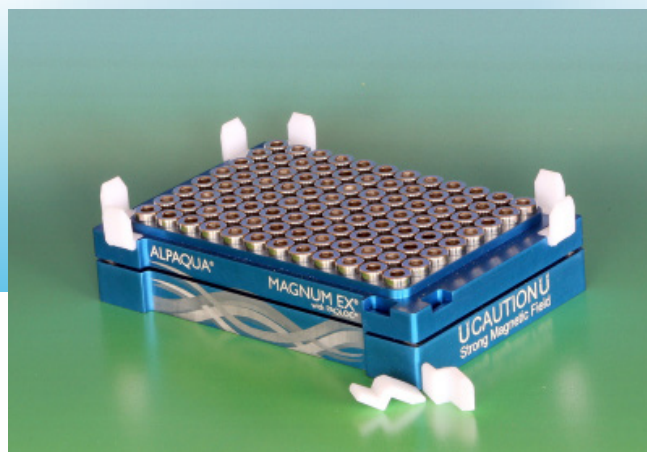


Use with springs (left) on the robot, or without springs on the bench



MAGNUM™ EX

2. Generation

Universal Magnet Plate

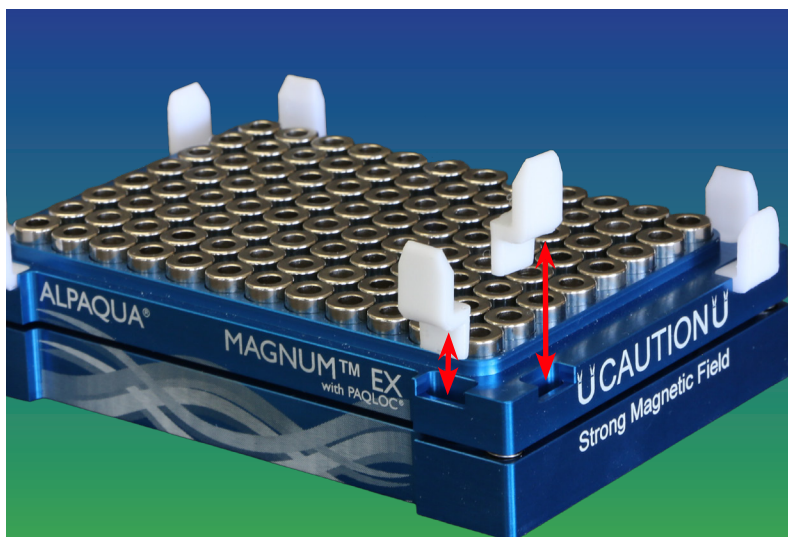

Eight corner tabs easily slide in and out of the plate.

There when you need them, out of the way when you don't.



PAQLOC Spring Lock:
Enable or disable the
spring cushion in seconds
with these four screws

Ordering Information

Buy online at www.alpaqua.com/shop, or place your order by phone, fax, or email at sales@alpaqua.com.



Product Name	Description	Catalog #
MAGNUM™ EX	Universal Magnet Plate 2. Generation (96-well)	A000580

Tel: (800) 690-1620

Outside US: +1 (978) 878-9489

Fax: (978) 865-9499

www.alpaqua.com

Alpaqua, PAQLOC, the wave logo, and Accelerating Genomic Discovery are registered trademarks of Alpaqua Engineering, LLC in the U.S. and/or other countries. Magnum is a trademark of Alpaqua Engineering, LLC.

*The PAQLOC Spring Cushion Lock is protected by U.S. Patent 12,372,497.

 **ALPAQUA®**
Accelerating Genomic Discovery®

©2025 Alpaqua Engineering. All Rights Reserved.
Document #BR000580-0925